

Accessories and Options for EP-TRACER



Measure the Difference.



SMART.
PRECISE.
RELIABLE.

SCHWARZER CARDIOTEK has a long history in the development of diagnostics for cardiology.


Today, our systems reflect the continuous cooperation with cardiologists and scientists, combining state-of-the-art research and development to deliver the smartest, most reliable and precise hardware and software solutions to our customers.

Measure the Difference.




Accessories and Options EP-TRACER

Catheter Connection Blocks and Accessories	4
ECG Cables	6
Pressure Cable Options for EP-TRACER	7
Suitcase	10
Keyboards and Mouse	11
Installation Components - Monitor Connection	12
Installation Components - Amplifier Connection	13
USB Connections, Powerbox and Filter	14
Order Form	15

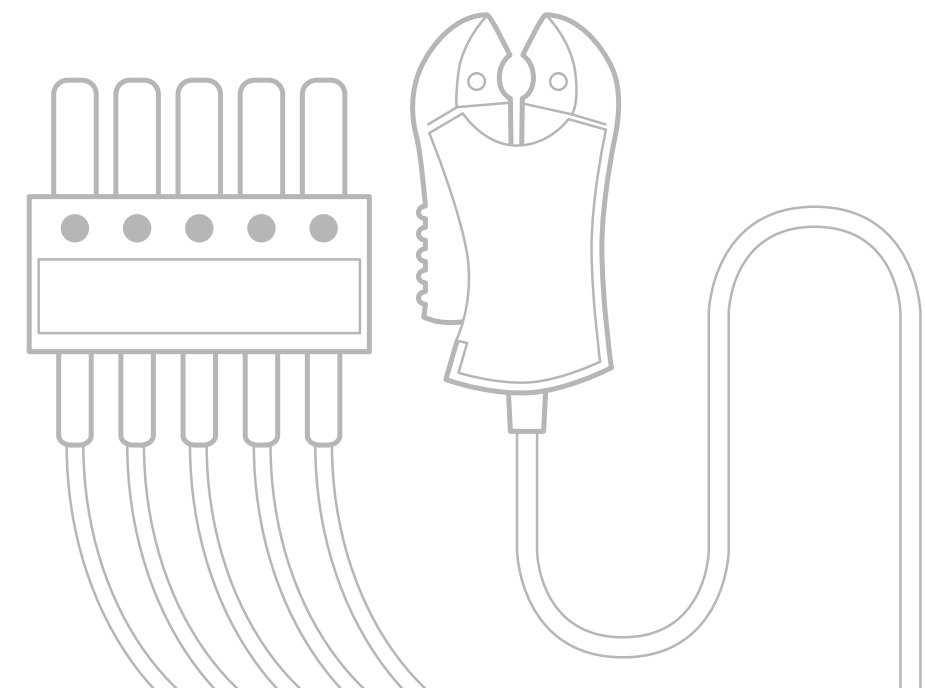

Catheter Connection Blocks

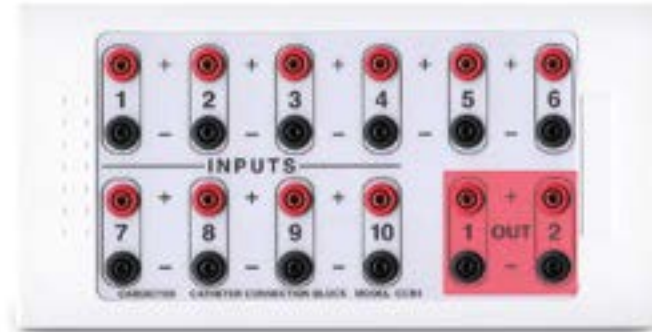


**01CCB04-G
CCB4G**

Catheter connection block with 20 intracardiac inputs for a maximum of 10 bipolar or unipolar channels. To be connected only in INT1 or INT2.




**01CCB05-G
CCB5G**

Catheter connection block with 16 intracardiac inputs for a maximum of 16 bipolar or unipolar channels. To be connected only from INT3 to INT6 (EP-TRACER 70 and 102 devices only).

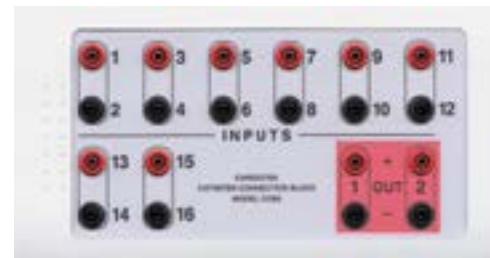

**01CCB06-G
CCB6G**

Catheter connection block with 10 intracardiac inputs for a maximum of 10 bipolar or unipolar channels. Specially designed for circular mapping catheters. To be connected only in INT1 or INT2.

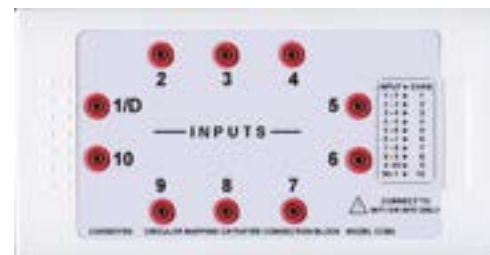

CCB Accessories



**03UBC02
UNI/BIP
Adapter Cables Set**

Splitter cables (set of 10 pieces). Used to split an intracardiac input to two different channels.



**03CCB01
CCB Cable**

Cable for connecting the catheter connection blocks to the EP-TRACER device. Suitable for all CCBs.



**06HLD01
Stainless Steel Holder
for 2 CCBs**

Stainless steel holder for two catheter connection blocks. Holder includes a clamp so it can be mounted to the norm rails of the patient bed.



ECG CABLES



Detail View



03ECG01 03ECG02 ECG Cable

10-electrodes cable to display 12-leads ECG.

03ECG01 Length 300 cm

03ECG02 Length 500 cm



Combine with



03ECG03 Translucent ECG Cable Set

Radiotranslucent extension cables for the ECG cable.

Specifications 03ECG03

Set with 7 pieces

Length 90 cm



Edwards Pressure Cable and Transducer



Length
390 cm

Detail View



03PRS01 Edwards Pressure Cable

Pressure cable for use with Edwards Lifesciences PX600F transducers.



Combine with



04PRS01 Edwards Pressure Transducer

Pressure transducer from Edwards Lifesciences model PX600F. Can be used only in combination with the Edwards pressure cable.

Specifications 04PRS01

Single-patient use

Sterile



Ohmeda Becton Dickinson Pressure Cable



Length
175 cm

Detail View



03PRS02 Pressure Cable for BD

Pressure cable for use with BD DTX plus transducers.



Combine with



04PRS02 Pressure Sensor BD DTX plus

Pressure transducer from BD model DTX plus. Can be used only in combination with the pressure cable for BD transducers.



Specifications 04PRS02

Single-patient use

Sterile

BBraun Pressure Cable



Length
295 cm

Detail View



03PRS04 Pressure Cable for BBraun

Pressure cable for use with BBraun transducers.



Combine with



BBraun Combitrans disposable transducer



Aluminium Suitcase for EP-TRACER 38



Detail View



08ACC01 Aluminium Suitcase for EP-TRACER 38

Aluminum suitcase. Used to safely store and transport an EP-TRACER 38 device plus accessories.



Length	Width	Height	Weight netto
500 cm	300 cm	15 cm	2,71 kg

Smart Keyboard



32KPD00 Smart Keyboard EP-TRACER

Smart Keyboard for single-key access to specific software functionalities.



32USK02 32USK03 Disinfectable Desktop Keyboard (D and US)

Types

32USK02 German version QWERTZ USB keyboard

32USK03 US version QWERTY USB keyboard



32OMS02 Disinfectable Desktop Mouse

Disinfectable mouse for EP-TRACER.



Monitor Connection



03DVI00 DVI Extender Set

Pair of DVI connectors to extend digital video signal by means of an optical fiber cable. The connectors have to be powered. Power cables are included.



03DPE00 Display Port Extender Set

Display Port extender connectors to extend digital video signal by means of an optical fiber cable. The connectors have to be powered. Power cables are included.



03FIB10 03FIB20 03FIB30 Optical Fiber

Can be used in combination with the DVI extender set, with the Display Port extender set or with the USB extender hub (page 13).

Types

03FIB10 Length 10.0 m

03FIB20 Length 20.0 m

03FIB30 Length 30.0 m



Amplifier Connection



03UEH00 USB Extender Hub

Pair of transmitter and receiver modules to extend USB communication by means of an optical fiber cable. Cables for powering the hubs are included.

To use with **03FIB10 / 03FIB20 / 03FIB30** (see page 12).

Combine with



03UEC20 03UEC30 Optical USB Extension Cable

Cables used to bridge larger distances for USB communication. Optical extender cables with fixed USB connectors.

Types

03UEC20 Length 20.0 m

03UEC30 Length 30.0 m



USB Connections, Powerbox & Filter



03USB01 USB Cable

USB cable of 1.8 meters to establish communication between the EP-TRACER device and PC. Also used with ProCart.

Length
180 cm



03USC01 USB Serial Converter

USB-to-serial converter. Used to connect serial or null-modem cables to any USB port of the EP-TRACER system.

Length
60 cm



03SNC01 Serial Cable for Cryo

Null-modem cable of 3 meters. Used to receive cryo ablation data from the Cryo Console from Medtronic.

Length
300 cm



05POB01 Powerbox

Power supply for the EP-TRACER device.



08FILO1 Filter Set

Set of 5 filters for fan (only for EP-TRACER 70 and 102 models).



ORDER FORM

Email customerservice@schwarzercardiotek.com

Phone +49 (0)7131 2774 500

Fax +49 (0)7131 2774 590

Your Contact Data

Hospital

Delivery Address

Department

Phone Number

Contact Name

Email

Order

Pricing Information

Page	Article-No. Schwarzer Cardiotek	Description	Order Quantity
CATHETER CONNECTION BOXES AND ACCESSORIES			
4	01CCB04-G	Catheter Connection Block CCB4-G	
4	01CCB05-G	Catheter Connection Block CCB5-G	
4	01CCB06-G	Catheter Connection Block CCB6-G	
5	03UBC02	UNI/BIP Adapter Cables Set	
5	03CCB01	CCB Cable	
5	06HLD01	Stainless Steel Holder for 2 CCBs	
ECG CABLES			
6	03ECG01	ECG Cable 3.0 m	
6	03ECG02	ECG Cable 5.0 m	
6	03ECG03	Translucent ECG Cable Set	
PRESSURE CABLE OPTIONS FOR EP-TRACER			
7	03PRS01	Edwards Pressure Cable	
7	04PRS01	Edwards Pressure Transducer	
8	03PRS02	Ohmeda Becton Dickinson (BD) Pressure Cable	
8	04PRS02	Ohmeda Becton Dickinson (BD) DTX plus Pressure Sensor	
9	03PRS04	BBraun Pressure Cable	
SUITCASE, KEYBOARDS AND MOUSE			
10	08ACC01	Aluminium Suitcase for EP-TRACER 38	
11	32KPD00	Smart Keyboard EP-TRACER	
11	32USK02	Disinfectable Desktop Keyboard D	
11	32USK03	Disinfectable Desktop Keyboard US	
11	32OMS02	Disinfectable Desktop Mouse	
INSTALLATION COMPONENTS MONITOR CONNECTION			
12	03DVI00	DVI Extender Set	
12	03DPE00	Display Port Extender Set	
12	03FIB10	Optical Fiber 10.0 m	
12	03FIB20	Optical Fiber 20.0 m	
12	03FIB30	Optical Fiber 30.0 m	
INSTALLATION COMPONENTS AMPLIFIER CONNECTION			
13	03UHE00	USB Extender Hub	
13	03UEC20	Optical USB Extension Cable 20.0 m	
13	03UEC30	Optical USB Extension Cable 30.0 m	
USB CONNECTIONS, POWERBOX AND FILTER			
14	03USB01	USB Cable	
14	03USC01	USB Serial Converter	
14	03SNC01	Serial Cable for Cryo	
14	05POB01	Powerbox	
14	08FIL01	Filter Set	

Date, Place

Signature

Hotline: +49 (0)7131 27 74 555
www.schwarzercardiotek.com

Schwarzer Cardiotek GmbH
Im Zukunftspark 3
74076 Heilbronn
Germany
Phone +49(0)7131 27 74 500
Fax +49(0)7131 27 74 590
info@schwarzercardiotek.com

CardioTek B.V.
Dr. Nolenslaan 107B
6136 GL Sittard
The Netherlands
Phone 1 +31 (0)85 48 67 64 0
Phone 2 +31(0)43 36 56 00 6